



✓ DESCRIPTION

EliGene® Coeliac RT is intended for the screening of alleles:

- HLA-DQ2 (DQA1*05/DQB1*02),
 - HLA-DQ8 (DQA1*03/DQB1*03:02),
 - HLA-DR4 (DRB1*04)
- from isolated DNA.

Specimens: Buccal swabs, blood
 Detection technology: TaqMan LNA

✓ PERFORMANCE STUDY

EliGene® Coeliac RT kit specificity was tested on 50 samples of human DNA with genotypes determined by DNA sequencing. Totally 50 samples were right determined by EliGene® Coeliac RT.

The clinical specificity of EliGene® Coeliac RT kit is 100%.

Analytical sensitivity: from 1 ng DNA molecules.

Specificity: 100 %


✓ ORDER DETAILS

Cat. No.: 90048-RT

50 reactions

✓ ADVANTAGES

- Detection of HLA alleles associated with Coeliac Disease
- Real-Time PCR
- Results in 60 min.
- Including internal and positive control

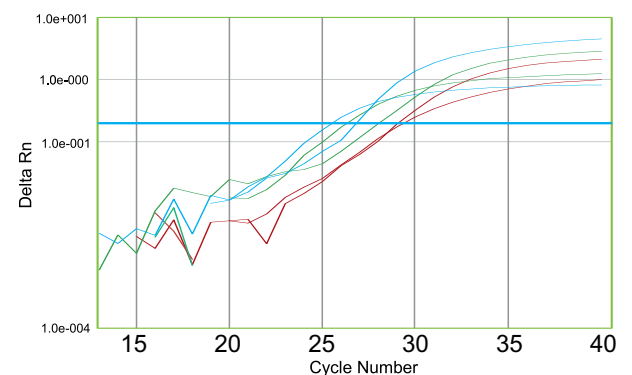

✓ COMPATIBLE DEVICES

ABI 7500 FAST (Applied Biosystems); LightCycler® 480 (Roche), RotorGene Q (Qiagen); Quant Studio 3 and 5 (ThermoFisher Scientific); CFX96 Touch Real -Time PCR Detection System (Bio -Rad)

CONTACT

ELISABETH PHARMACON
 ROKYCANOVA 4437/5, BRNO-ŽIDENICE
 CZECH REPUBLIC | +420 542 213 851
www.elisabeth.cz

